


SEDEX Low Temperature Evaporative Light Scattering Detector

SAMPLE: CORTICOSTERONE

Micro HPLC

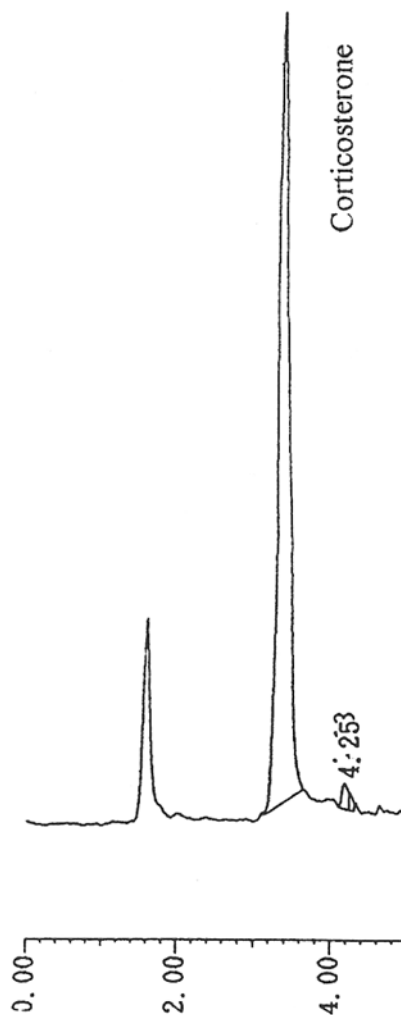

Limit of detection: 10 ng

CONDITIONS

Column	<i>Spherisorb S5W</i> 250 x 2 mm
Injected Volume	200 ng/ μ l
Eluent	Chloroform/Water (96/4)
Flow Rate	150 μ l/min
Pressure	2.4 bars
Gas	Air
Temperature	70° C

Analyzed By : M. Kawakatsu and H. Kotagoniguichi, MS Trading, Osaka JAPAN