


**UNDERIVATIZED AMINO ACIDS**

**HPLC**

**WITH EFFECT OF ION PAIRING AGENT**

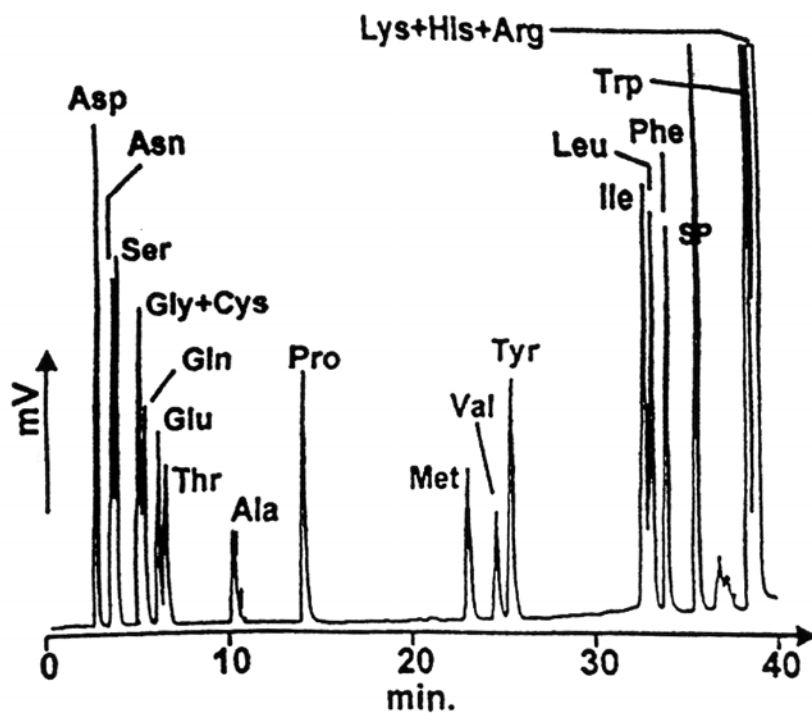

Silica-C8 ( 125X4.6 mm ID)  
 elution gradient: H<sub>2</sub>O (+volatile ion pairing agent)  
 0% to 40% organic solvent in 40 min.  
 flow rate : 1 ml/min

**CONDITIONS**

Pressure	2 bars
Gas	Air
Gain	10
Temperature	45°C

Analyzed by : Mr Dreux et al. - UNPUBLISHED RESULTS - SP : system peak